

# Uni-Shunt® System


The Uni-Shunt® System with Reservoir incorporates a double dome access port to facilitate injections and aspirations of CSF samples. It is a continuous length of barium-impregnated silicone tubing with an access reservoir made of self-sealing silicone which can be punctured with a 25 gauge or smaller Huber type needle.

The simplicity of the one-piece design means that there are no connectors, so the possibility of disconnection and subsequent complications is eliminated.

The catheter has X-ray detectable tantalum markings and tantalum tips on the ventricular and peritoneal ends. X-ray detectable pressure markings at the distal tip indicate the slit valve closing pressure.

The Uni-Shunt Catheter with Reservoir Kit is supplied sterile with the following:

- Blunt tip needle, 17 gauge, disposable, for flushing and testing purposes
- Disposable malleable catheter passer, short. Each 90 cm and 102 cm length kits also include one disposable malleable catheter passer, long
- Uni-Shunt Catheter with reservoir anchoring clips, 2
- Uni-Shunt Catheter with reservoir introducer

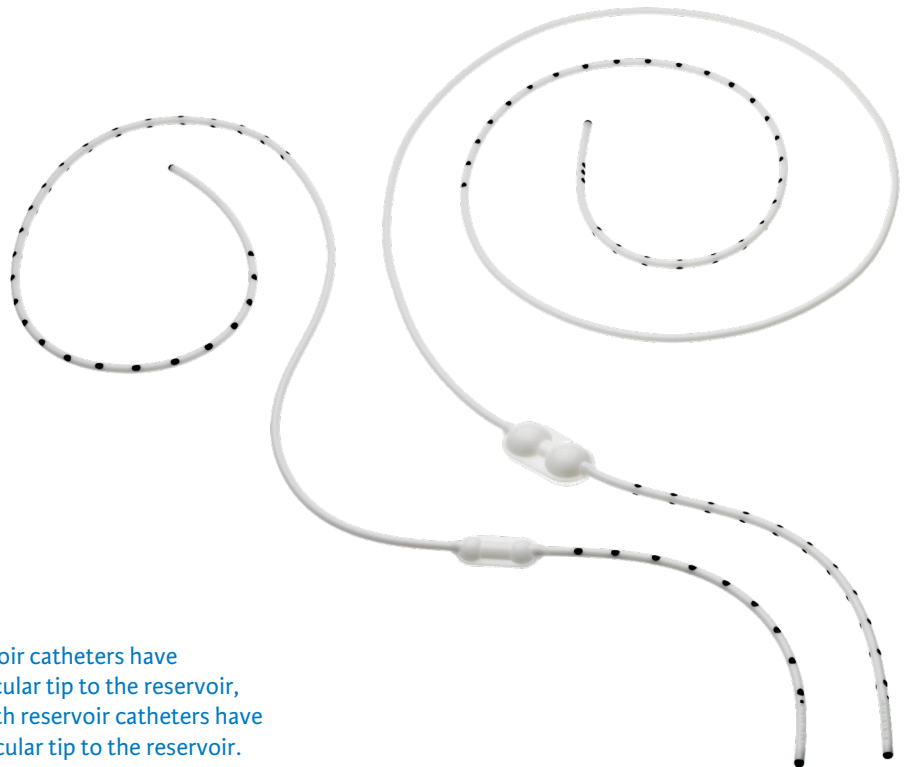

**Note:**

The 53 cm and 62 cm length reservoir catheters have a distance of 11 cm from the ventricular tip to the reservoir, the 78 cm, 90 cm, and 102 cm length reservoir catheters have a distance of 13 cm from the ventricular tip to the reservoir.



Pediatric Reservoir Kit

**Uni-Shunt System**

Catalog#	Length	Closing Pressure (Measured in cm of water)
82-8501	53 cm	2-5 cm H <sub>2</sub> O Low
82-8502	53 cm	5-9 cm H <sub>2</sub> O Medium
82-8511	62 cm	2-5 cm H <sub>2</sub> O Low
82-8512	62 cm	5-9 cm H <sub>2</sub> O Medium
82-8521	78 cm	2-5 cm H <sub>2</sub> O Low
82-8522	78 cm	5-9 cm H <sub>2</sub> O Medium




Adult Reservoir Kit

Catalog#	Length	Closing Pressure (Measured in cm of water)
82-8531	90 cm	2-5 cm H <sub>2</sub> O Low
82-8532	90 cm	5-9 cm H <sub>2</sub> O Medium
82-8533	90 cm	9-14 cm H <sub>2</sub> O High
82-8541	102 cm	2-5 cm H <sub>2</sub> O Low
82-8542	102 cm	5-9 cm H <sub>2</sub> O Medium
82-8543	102 cm	9-14 cm H <sub>2</sub> O High

**Holter® Lumbo-Peritoneal Catheters**

Manufactured as a one-piece silicone rubber with an ID of 1.2 mm and an OD of 2.5 mm. The 3 cm T-portion has four sets of three proximal holes, and has 1.0 mm ID and 1.8 mm OD. The peritoneal catheter end has three longitudinal slits (120° apart) to regulate flow.

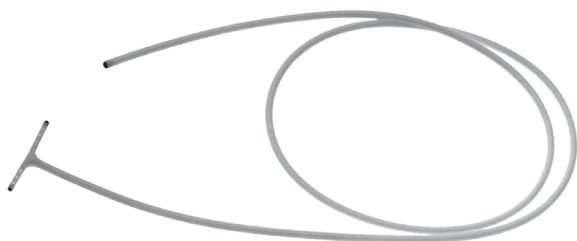

"T" Detail

Catalog#	Description
82-1698	Holter® Lumbo-Peritoneal Catheter, T-shaped, one-piece <ul style="list-style-type: none"> <li>• Length 65 cm</li> <li>• 3.5 mm slit valves</li> </ul>

Availability of these products might vary from a given country or region to another, as a result of specific local regulatory approval or clearance requirements for sale in such country or region.

- Non contractual document. The manufacturer reserves the right, without prior notice, to modify the products in order to improve their quality.
- Warning: Applicable laws restrict these products to sale by or on the order of a physician.
- Consult product labels and inserts for any indication, contraindications, hazards, warnings, precautions, and instructions for use.

---

**For more information or to place an order, please contact:**

**United States, Canada, Asia, Pacific, Latin America**

USA 800-654-2873 ■ 888-980-7742 fax

International +1 609-936-5400 ■ +1 609-750-4259 fax

[integralife.com](http://integralife.com)

