



Limit uncertainty with a thinner matrix to meet complex needs.

Integra® Wound Matrix (Thin) combines the proven technology of our Integra Wound Matrix with the versatility of a thinner matrix. This product can be applied to complex anatomical sites that require a more conformable graft.

Engineered for Wounds

- Integra Wound Matrix (Thin) is approximately 0.4mm thick, half the thickness of Integra Wound Matrix
- Supports a healing environment for wound management, including partial and full-thickness wounds and donor sites
- Cross linked collagen-chondroitin-6-sulfate

Proven Performance


- Over 30 years of science and innovation in the development of collagen technology
- Integra's collagen is the base material of implants used successfully in over 12 million procedures

Description

Integra Wound Matrix (Thin) is a collagen-glycosaminoglycan wound dressing that maintains and supports a healing environment for wound management. Integra Wound Matrix (Thin) is supplied sterile and is intended for one-time use.

Benefits

- Provides immediate wound coverage
- Highly conformable for various anatomical sites
- Room temperature storage and 24 month shelf life
- Terminally sterilized
- Available in sizes specifically designed for partial and full-thickness wounds and donor sites

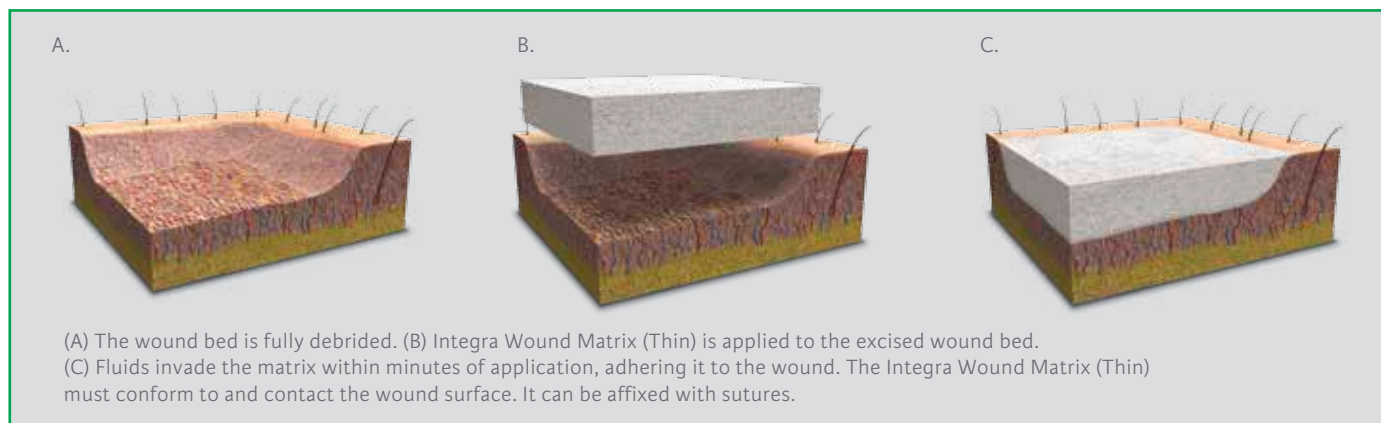

Indications

Integra Wound Matrix (Thin) is indicated for the management of wounds including: partial and full-thickness wounds, second degree burns, pressure ulcers, venous ulcers, diabetic ulcers, chronic vascular ulcers, tunneled/undermined wounds, surgical wounds (donor sites/grfts, post-Moh's surgery, post-laser surgery, podiatric, wound dehiscence), trauma wounds (abrasions, lacerations, skin tears) and draining wounds. The device is intended for one-time use.

Contraindications

This device should not be used in patients with known sensitivity to bovine collagen or chondroitin materials. The device is not indicated for use in third degree burns.

Precautions

The following complications are possible with the use of wound dressings. If any of the conditions occur, the device should be removed: infection, chronic inflammation (initial application of wound dressings may be associated with transient, mild, localized inflammation), allergic reaction, excessive redness, pain or swelling.

Integra Wound Matrix (Thin)

Catalog Number	Size	Units/Case
52021T	2in x 2in (5cm x 5cm)	1 Sheet/Box
54051T	4in x 5in (10cm x 12.5cm)	1 Sheet/Box
54101T	4in x 10in (10cm x 25cm)	1 Sheet/Box

See package insert for full prescribing information.

For more information or to place an order, please contact:
 Integra ■ 311 Enterprise Drive, Plainsboro, NJ 08536
 800-997-4868 USA ■ 609-936-5400 outside USA ■ 888-980-7742 fax
integralife.com

